

SOLID PLUS

WWW.ROTAID.COM



 **ROTAID™**
AED CABINETS

ROTAID SOLID PLUS

PRODUCT SHEET

The Rotaid Solid Plus is an innovative AED wall cabinet that ensures a recognizable and robust storage of your AED device. The characteristic round shape and green colour, make it easy to be found and recognized and comes with a high quality appearance. The large AED logo, positioned in the transparent centre of the cover, leaves no room for interpretation. This is a life saving device.

The cabinet can be opened by an intuitive twist that removes the whole cover. The responder therefore has full access to the AED that hangs via a bracket for a quick and easy grab. The cover is attached to the back plate via a belt and needs no further attention from the responder. A stainless steel hook, attached to the bottom of the cabinet, can hold a seal, tag or other appliance.

The cover, equipped with recesses for an optimal grip and directive arrows, is made from 100% Polycarbonate. With a thickness of 5 millimetres it offers a solid protection and ensures a high grade of water- and dust resistance. Furthermore, the cover is fitted with UV filters, to protect from direct and indirect sunrays. Appearance decay, rust and accumulation of residue are excluded.



THE ROBUST DESIGN AND THE HIGH DEGREE OF DUST AND WATER DENSITY OFFER MAXIMUM PROTECTION.



100% DURABLE, BECAUSE OF MODERN PRODUCTION TECHNIQUES AND THE USE OF RECYCLED PLASTICS




FULL ACCESS TO THE AED


LEATHER BELT HOLDS COVER WHEN OPENED

ROTAID SOLID PLUS

TECHNICAL DETAILS

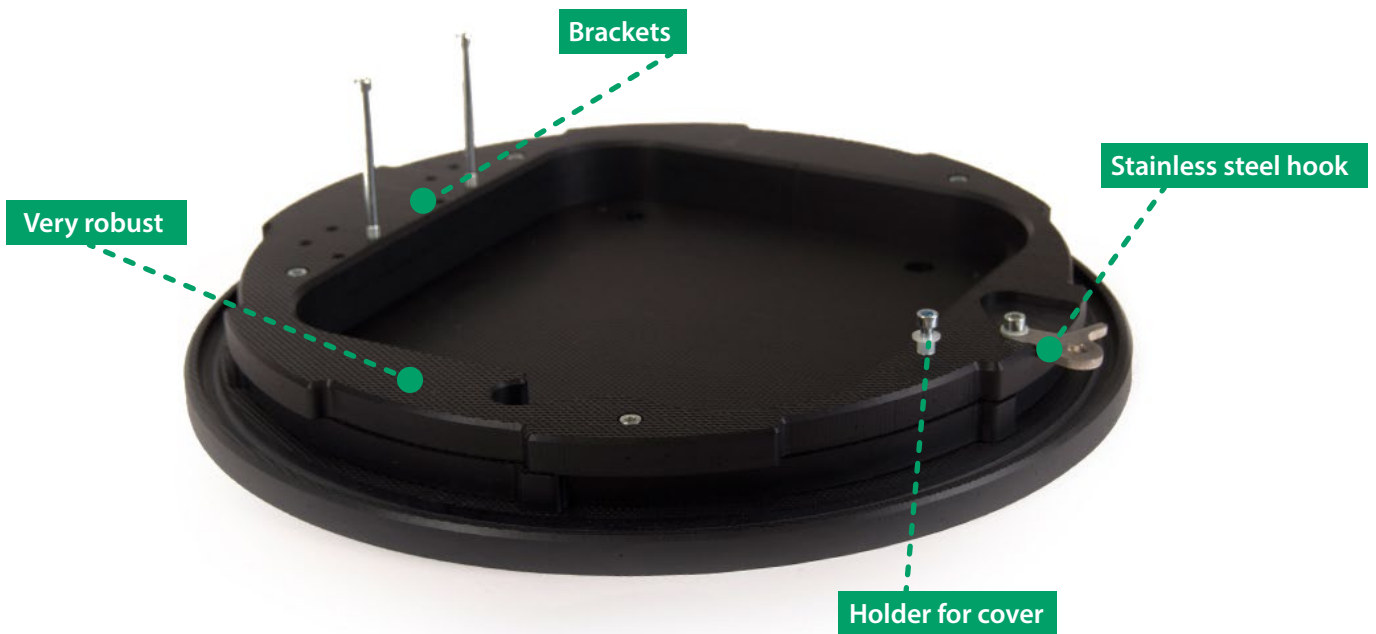

The back plate is made from 100% HDPE plastics and offers space and protection to practically all types of AED devices on the market.

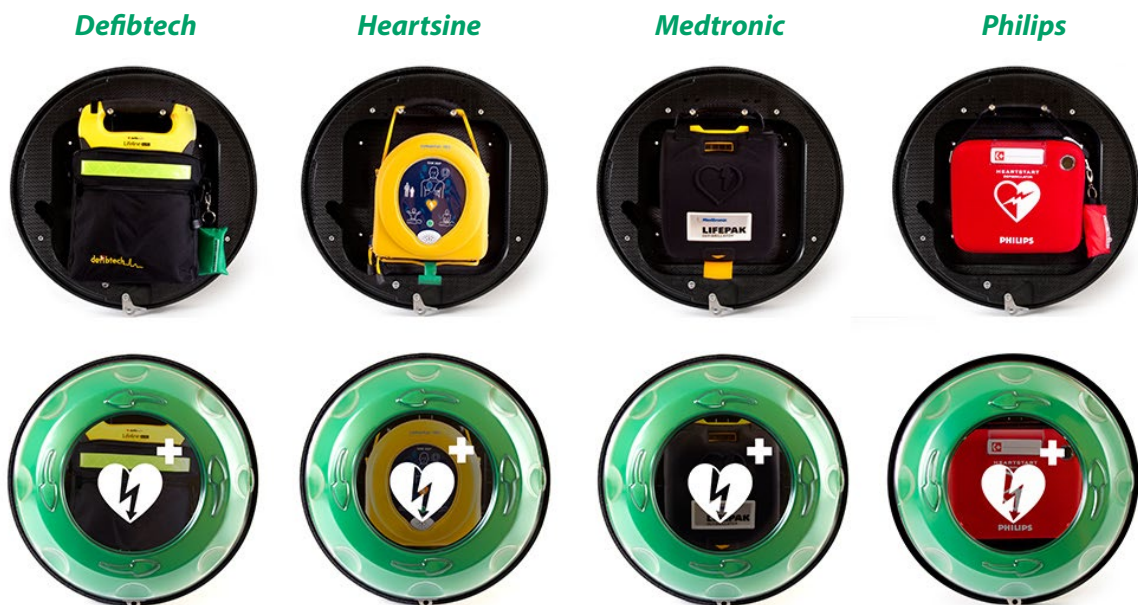

Also compatible with *Cardiac Science, CU Medical, Primedic, Schiller and Zoll.*

Dimensions	Length 48 cm Width 48 cm Depth 18 cm Weight 6,7 kg
	Length 18,9" Width 18,9" Depth 7,1" Weight 14,7 lbs

Packed individually by box, fastening materials included.